# Gebruiksaanleitung Manual Mode d'emploi

(D)

Seite 2

(GB)

Page 5

(F)

Page 8

**Alpha<sup>®</sup> Aufrichthilfe Standard / Magic**  
**Alpha<sup>®</sup> Raising Aid Standard / Magic**  
**Alpha<sup>®</sup> Sangle Lève-malade Standard / Magic**

Magic

Standard

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## 1. Specifications and Scope of Supply

**Item no.:** PM-3015 Magic  
**Dimensions:** 65 x 25 cm  
**Material:** Nylon, latex free  
**Color:** blue/red  
**ILN:** 4250355330154  
**Storage:** -5 – 40° C  
 protect from sunlight  
**Patent no.:** EP 2248502

### Scope of supply:

- Raising aid Magic with flexible safety harness and removable core
- Manual

**Item no.:** PM-3012 Standard  
**Dimensions:** 60 x 20 cm  
**Material:** Nylon, latex free  
**Color:** blue/red  
**ILN:** 4250355330123  
**Storage:** -5 – 40° C  
 protect from sunlight



### Scope of supply:

- Raising aid Standard with removable core
- Manual

### Accessory:

PM-3010 safety loop belt for supporting the patient transfer of persons with limited body tension / limited space available

Made in Germany
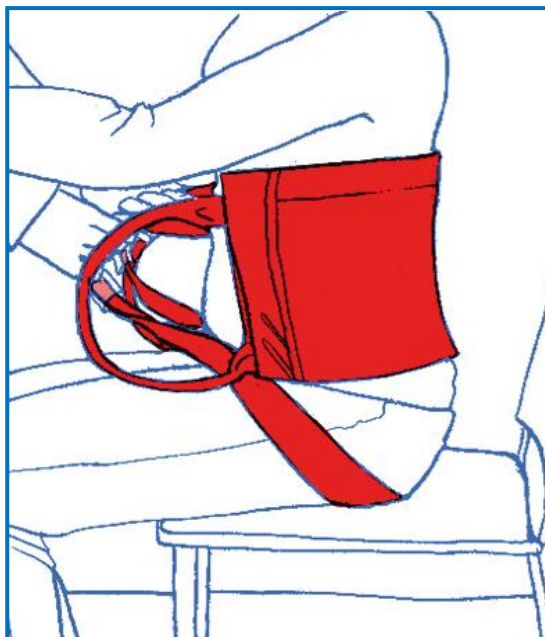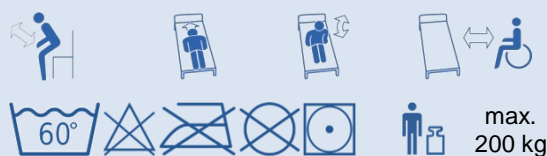
## 2. Use / Placement

The raising aids are universal aid with many potential applications.


### 2.1 Raising Aid Magic

The raising aid Magic is used to move patients within the bed or transfer them from a bed to a wheelchair. It is ideal for use in conjunction with the sliding board when transferring the patient, e. g. to a wheelchair.

You can use the raising aid Magic with safety harness on top or downward depending on the situation. The additional flexible harness finds its position by moving it repeatedly to the right and to the left behind the patient's back in sitting position. It is also an option to place the belt behind the patient with both hands under the patient's bottom. When using the additional harness, please use the handle that is connected to the harness.

## Example of use:
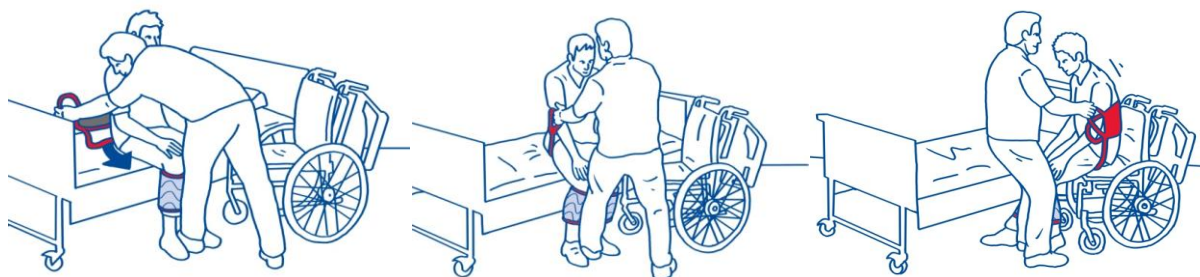
### 2.2 Raising Aid Standard

The raising aid Standard is used to move patients within the bed or transfer them from a bed to a wheelchair. It is ideal for use in conjunction with the sliding board when transferring the patient, e. g. to a wheelchair.

### 2.3 Accessory (option): safety loop belt for patient support

The safety loop belt can be used optionally for patients with limited body tension.

The safety loop belt enables also especially the safe transfer for patients that sit in very narrow wheel chairs or seat pans that prevent applying the flexible holding belt of the raising aids under the patient. Both loops of the safety loop belt have to be attached to the grips of the holding belt and with help of the clip locks then applied under the patient's legs.

The safety loop belt is exclusively designed for the short-term use during patient transfers and only in combination with the PM-Alpha products (holding belt and raising aid Magic / Standard)

**These products are suitable for re-use.**

## 3. Maintenance / Care

---

The product does not require any special maintenance. It can be washed in a washing machine at 60° C and is suitable for low dryer temperatures (45 – 50° C). When using non-standard washing agents, bleaching agents, softener or other chemical washing additives, the functional properties of the product can strongly degrade. The robust material can be disinfected with the usual wipe disinfectants. Check the compatibility before use where required on the basis of the disinfectant's instructions. We recommend the use of alcohol- and aldehyde-free agents, e.g. *Trionic®*. Before washing, remove the core of the product and wipe it or disinfect it if necessary. There is also some care information on the label.

## 4. Warranty

---

- Petermann gives a warranty that the product is devoid of material defects and manufacturing defects.
- The warranty is valid for a period of **3 years**. It does not cover any defects which are due to improper use.
- Petermann cannot be held liable for any accidental or indirect injury or damage.

## 5. Visual and Function Inspection / Safety Information / Re-use

---

Repositioning products are also subject to natural aging. Check the equipment regularly and carefully to make sure it is working properly and especially before re-use. Make sure that the product is in good condition and that the material and the seams are not damaged. Check the function of the product by applying pressure and loading to the handles and make sure that the material resists the load. The high-friction has a rough surface to prevent slippage. The best way to test these properties is by comparing the product with a new product.

For re-using the product it must be washed and disinfected before, see 3.

### **WARNING:**

- ***Faulty products must no longer be used.***
- ***If the patient's mobility is restricted, it is essential to supervise the motions of transfer by providing trusted help.***
- ***When using non-standard washing agents, bleaching agents, softener or other chemical washing additives, the functional properties of the product can strongly degrade.***
- ***Carefully check the transfer surfaces, making sure they are stable and brakes are applied (e. g. wheelchair, bed).***
- ***Use the safety loop belt exclusively for short-term patient transfer and only in combination with Alpha holding belt and raising aid Magic/Standard.***
- ***Mind the possible hazards of tripping and getting caught by the belt systems!***

## 6. Product Labeling

---

The Petermann Alpha<sup>®</sup> raising aids are medical products, class I and complies with the MDD directive 93/42/EEC.

Petermann GmbH is certified in accordance with the directive DIN EN ISO 13485 by the German TUV.

## 7. Disposal Information

---

For the raising aids there are no special disposal rules; they can be disposed of in accordance with national and local regulations as necessary.

If you have any questions about the use of our products, please call our service number:

**+49 (0) 9868 9339-0**

Petermann GmbH • Ziegelhaus 12 • D-91601 Dombühl  
Postal address: Petermann GmbH • Waldweg 3 • D-91601 Dombühl  
Phone: +49 (0) 9868 9339-0 • Fax: +49 (0) 9868 9339-20  
E-mail: petermann@pm-med.de • www.pm-med.de

"Information subject to technical modifications to our products and deviations in colour or printing errors"

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Petermann GmbH est certifié par le TUV conforme à la directive DIN EN ISO 13485.
Adresse postale: Petermann GmbH • Waldweg 3 • D-91601 Dombühl